

# mediTECH

MEDICAL GAS EQUIPMENT

## 4 Way Flow Selector On Stand Mount

BS Schrader Probe on 2M Supply Hose



- **No Calibration required EVER**
- **Lightweight construction**
- **Available with International Fittings**
- **7 Year Warranty**
- **7 Year service Interval**
- **Mass casualty device**

Excellent performance and robustness, coupled with a very lightweight design. Manufactured from a strong aluminium stand mount. Supplied with 2M supply hose to connect 4 BAR Oxygen source (longer hose length available upon Request).

For use in all Hospital and Emergency Environments, the 4 Way Flow Selector is designed to be plugged into Medical Oxygen Wall Outlets.

The Meditech 4 Way Flow Selector has been designed to deliver precision flow rates, without the need for any calibration. It has a standard 2M hose length (other hose lengths available upon request.)

Supplied with Push Fit “fir-tree” for direct attachment of Oxygen Therapy accessories.

**AVAILABLE GAS TYPES:**  
OXYGEN | MEDICAL AIR | ENTONOX

### THERAPY OUTLET TYPES



FIR TREE



9/16<sup>th</sup> THREADED  
THERAPY CONNECTOR



3/8" BSP THREADED  
THERAPY CONNECTOR




1/4" BARB

### PROBE TYPE



BS PROBE

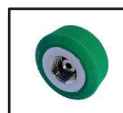

DISS



DIN



AGA




AFNOR



SIS



CZECH

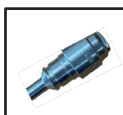

UNIFOR



# SPECIFICATION

Intended Gas: .....

Medical Oxygen (O<sub>2</sub>)

Inlet pressure: .....

4 bar (400 kPa)

Set Flow Rates in lpm: .....

Three flow rates options available:

1/2, 1, 2, 3, 4, 6, 8, 10, 15 lpm

1/2, 1, 2, 3, 4, 6, 8, 10, 15, 25 lpm

10, 20, 30, 35, 40, 50, 55, 60, 70 lpm

Inlet Connections: .....

Schrader (BS 5682:2015) and other international connections

Outlet (Patient) Connections. ....

Firtree style barbed outlet.

Intended Life: .....

10 Years

Service Interval: .....

7 Years

Warranty Period: .....

7 Years

Complies With: .....

ISO 15002: 2023

Weight including hose: .....

1.6Kg

Dimensions: .....

Height: 15CM x 12CM (excluding hose)

## TO PURCHASE OR FOR FURTHER DETAILS:

**Contact us by email or telephone for more details, including specific part numbers, prices, etc.**

B.N.O.S. Meditech Ltd.

9 Fifth Avenue, Bluebridge Industrial Estate,

Halstead, Essex CO9 2SZ, England, UK.

Tel: +44 (0)1787 479475 • Fax: +44 (0)1787 477747

E-mail: [sales@meditech.uk.com](mailto:sales@meditech.uk.com) • [www.meditech.uk.com](http://www.meditech.uk.com)

Meditech Flow Selectors are CE marked in accordance with the European Medical Devices Directive 93/42/EEC and UKCA marked accordance with the UK Medical Devices Regulations 2002. All items in this range are sold with a return to manufacturer warranty against manufacturing defects. These products are part of a series of medical gas and resuscitation products available from B.N.O.S. Meditech. If the product you require is not shown here, or if you require more information, please contact your distributor or contact us directly. Specifications may be subject to change without notice.

B.N.O.S. Meditech Ltd Terms and Conditions of sale apply, available on application or from our website.



# medITECH

MEDICAL GAS EQUIPMENT

Distributor: