

mediTECH
MEDICAL GAS EQUIPMENT

microVENT™ World

Paramedic Resuscitation Kit



- **Automatic & Manual resuscitation.**
- **Pneumatic Oxygen powered resuscitator.**
- **For adult & child use over 10kg.**
- **Hand held, compact & easy to use.**

The microVENT World Paramedic Resuscitation Kit is designed for professional emergency services when breathing has stopped, when breathing assistance is needed or if oxygen therapy is required.

At the heart of the kit is the Microvent World resuscitator. This is an oxygen powered resuscitator that can be set to automatic or used manually. The microvent has a remarkably efficient use of oxygen, as long as 45 minutes or more from one cylinder.

The case is reinforced and padded to hold and protect the equipment. It is fitting with twin carrying handles and a shoulder strap. The resuscitator can be swiftly carried to the casualty.

The oxygen regulator fits to pin index medical oxygen cylinders. It has a BS Schrader connection to power the microvent. There is also a multiple flow therapy outlet for supplying the therapy mask or to use with other equipment needing a constant oxygen supply.




microVENT™ World Resuscitator

Part No. 670-0531-00

A hand held, light weight resuscitator with compact dimensions, the microVENT™ World has ergonomic styling to give easy operation with controls in a natural position. It is constructed with a tough reinforced plastic casing to withstand use in difficult situations.

The microVENT™ World patient valve features 22/15mm connections to take standard airway devices and face masks. The oxygen supply hose is 2 metres long and made of resilient reinforced plastic to retain flexibility in cold conditions. A British Standard (BS) medical oxygen supply probe enables use with all BS standard portable and piped oxygen supplies. Other international hose fittings are also available.

The patient assist valve on advanced models synchronises the resuscitator whilst in automatic mode with patient demand. Airmix models supply 100% or approx. 50% O2 concentration.

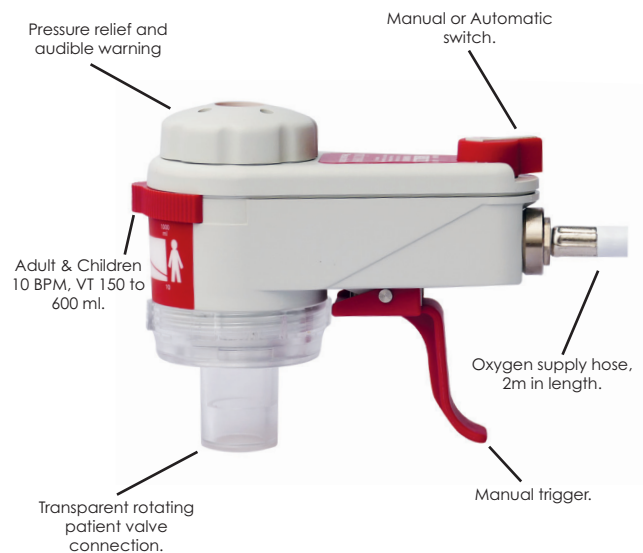

microVENT™ World Paramedic Resuscitation Kit

Part No. 670-0130-01



Please note that the Oxygen Cylinder is NOT included within the kit. This needs to be purchased separately. Please contact a member of our sales team for more info.

Medical Oxygen Pin Index Regulator

Part No. 670-0484-00



The microVENT™ World Resuscitation Kit includes a Meditech Oxygen Therapy Regulator with a Pin Index Inlet connector and a BS Medical Oxygen Schrader outlet. It is also fitted with a multiple therapy flow selector and therapy outlet.

The flow rates available are:

1/2, 1, 2, 3, 4, 6, 8, 10 and 15 litres per minute.

Accessories



Hook Ring
Part No. 673-0444-00



Head Harness
Part No. 673-0445-00



Manual SUCTION UNIT
Part No. 673-0448-20

TO PURCHASE OR FOR FURTHER DETAILS:

Contact us by email or telephone for more details, including specific part numbers, prices, etc.

B.N.O.S. Meditech Ltd.

9 Fifth Avenue, Bluebridge Industrial Estate,
Halstead, Essex CO9 2SZ, England, UK.

Tel: +44 (0)1787 479475 • Fax: +44 (0)1787 477747

E-mail: sales@meditech.uk.com • www.meditech.uk.com

CE
2797

UK
CA
0086

Components manufactured by BNOS Meditech and supplied within the kit are CE marked in accordance with the European Medical Devices Directive 93/42/EEC and UKCA marked in accordance with the UK Medical Devices Regulations 2002. All items in this range are sold with a return to manufacturer warranty against manufacturing defects. These products are part of a series of medical gas and resuscitation products available from B.N.O.S. Meditech. If the product you require is not shown here, or if you require more information, please contact your distributor or contact us directly. Specifications may be subject to change without notice.

B.N.O.S. Meditech Ltd Terms and Conditions of sale apply, available on application or from our website.

medITECH
MEDICAL GAS EQUIPMENT

Distributor: