


BEDSIDE MEDICAL OXYGEN SYSTEM

Supplied by BNOS Meditech Ltd

OXYGEN SIGN
674-0106-00

TWIN FLOW SELECTOR
0-15LPM = 710-0101-22
0-25lpm = 710-0101-23

**PIN INDEX BS SCHRADER
OUTLET**
673-0001-10

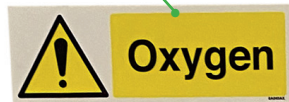

CYLINDER RACK
674-0107-00

HUMIDIFIER
Single use - 674-0047-00
X2 UNITS SUPPLIED PER KIT

OXYGEN THERAPY CONSUMABLES
High concentration face mask (adult)
673-9006-00
X10 UNITS SUPPLIED PER KIT

**40L MEDICAL OXYGEN
CYLINDER**

- Filled (676-0001-00)
- Empty (676-0002-00)

SERVICING REQUIREMENTS

Require 1 service at 5 years:

- 680-0001-10
- 680-0501-10
- 710-0101-23

Change after 3 years:

- 674-0047-00
- 673-9006-00



DESCRIPTIONS:

Meditech's Medical Oxygen Systems meet the latest IMO IMDG / MFAG regulations for shipboard use to supply Oxygen to two persons at the same time.

In the IMO Bedside solution, the Oxygen Cylinder is placed between two beds, mounted against the bulkhead with the cylinder rack. The regulator with a fixed flow of 50L/min is mounted on the top of the Cylinder and the Flow Selector is mounted on the bulkhead. The double Oxygen Outlet is fixed to the wall and the hose is connected to the Regulator. Two Flow Selectors, with quick connections, are placed in the double Oxygen Wall Outlet, which enables the Oxygen flow to be regulated between 0-15L/min. A humidifier is connected to each Flow Selector and this has a disposable Oxygen Therapy Mask attached to it.

- All equipment is supplied in one box ready for easy installation
- Installation can be completed without any trouble during a voyage and can easily be performed by the crew.
- All equipment is available through our extensive network.
- Future service is supported by our service organisation including replacement of components
- Compatible with Medical Oxygen Cylinders

Approvals:

Meets current IMO | IMDG | MFAG regulations

medITECH
MEDICAL GAS EQUIPMENT

TO PURCHASE OR FOR FURTHER DETAILS:

Contact us by email or telephone for more details, including specific part numbers, prices, etc.

B.N.O.S. Meditech Ltd.

**9 Fifth Avenue, Bluebridge Industrial Estate,
Halstead, Essex, CO9 2SZ, England, UK.**

TEL: +44 (0) 1787 479 475

FAX: +44 (0) 1787 477747

EMAIL: sales@meditech.uk.com

WEBSITE: www.meditech.uk.com

These products are part of a series of medical gas and resuscitation products available from B.N.O.S. Meditech. If the product you require is not shown here, or if you require more information, please contact your distributor or contact us directly.

Specifications may be subject to change without notice. B.N.O.S. Meditech Ltd Terms and Conditions of sale apply, available on application or from our website.

Distributor: